


## Adapter strand Package

with green dye **Ga**

Dynamic Biosensors GmbH & Inc.  
ASP-1-Ga v5.1

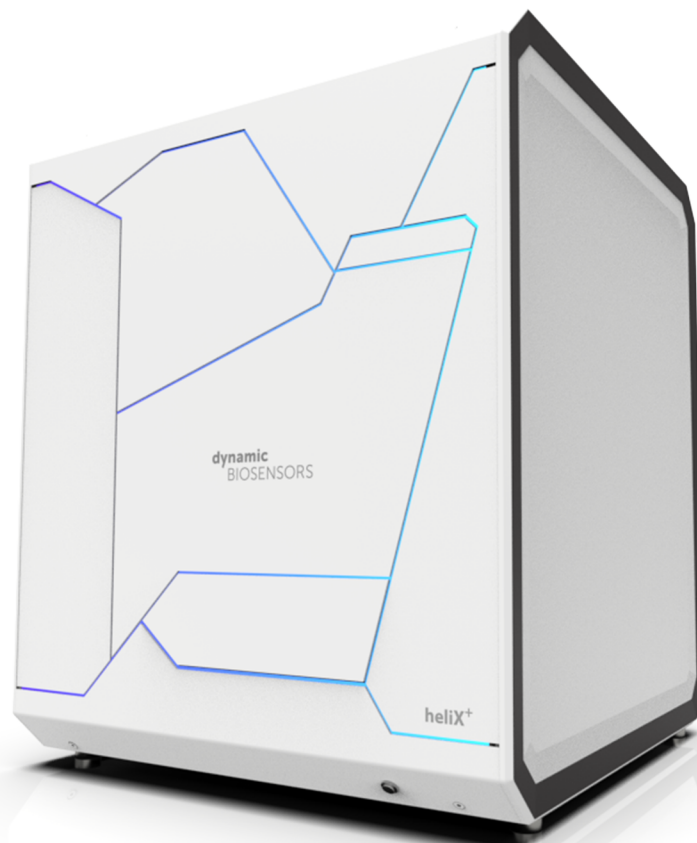


## Key Features

- **Adapter strand 1 - Ga** and **Adapter strand 2 - Ga - lfs** (ligand free strand) for functionalization of **heliX<sup>®</sup> Adapter Chip** on Spot 1 and Spot 2, respectively.
- Compatible with **heliX<sup>®</sup> Adapter Chip**.
- Includes **Adapter strands** for **50 regenerations**.
- Ideal for **MIX&RUN** sample preparation.
- **Adapter strands 1** and **2** carry a moderately hydrophilic green dye (**Ga**) with a single negative net charge.

## heliX<sup>®</sup> Adapter Chip Overview

2 spots with 2 different anchor sequences for DNA-encoded addressing.

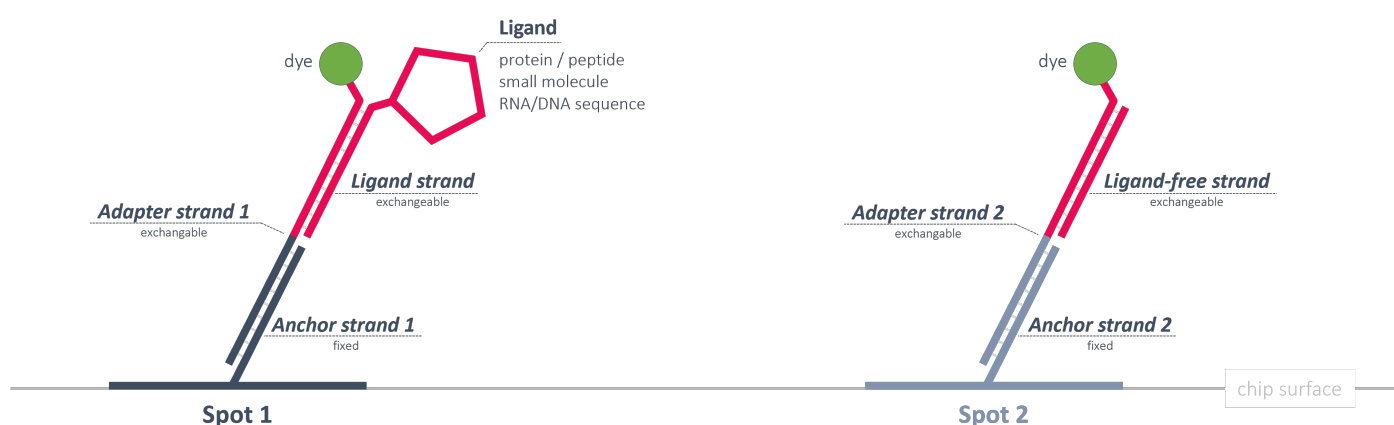

## Product Description

Order Number: **ASP-1-Ga**

Table 1. Contents and Storage Information

Material	Cap	Concentration	Amount	Buffer	Storage
<b>Adapter strand 1 - Ga</b>	Black	400 nM	5 x 100 µL	TE40 <sup>[1]</sup>	-20°C
<b>Adapter strand 2 - Ga - lfs</b>	White	200/250 nM	5 x 200 µL	TE40 <sup>[1]</sup>	-20°C

For research use only.

This product has a limited shelf life, please see expiry date on label.

To avoid many freeze thaw cycles please aliquot the nanolever.

## Preparation | MIX&RUN

In-solution hybridization of adapter and ligand strands:

1. Mix **Adapter strand 1 - Ga** (400 nM) and conjugated **Ligand strand** (500 nM) at 1:1 ratio (v/v).
2. Incubate the solution of step 1 at **RT** at **600 rpm** for **30 min** to ensure complete hybridization.
3. Mix solution of step 2 and **Adapter strand 2 - Ga - lfs** (200 nM) at 1:1 ratio (v/v).

Solution is ready to use for biochip functionalization.

Stability of the solution is related to the stability of the ligand molecule.

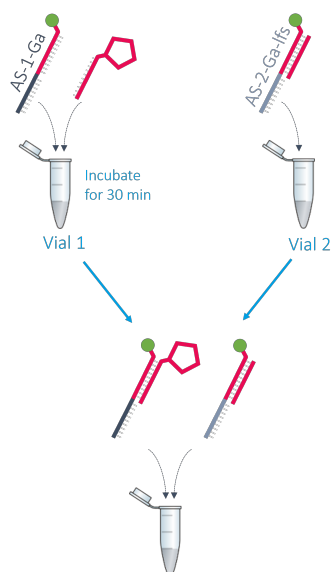


Table 2. Additional material for functionalization of spot 1 and reference spot 2.

Material	Concentration	Buffer	Related Product Name	Order No
<b>Ligand strand</b> carrying the conjugated ligand	500 nM	PE40 <sup>[2]</sup>	heliX <sup>®</sup> Amine Coupling Kit 1	HK-NHS-1

## Example

Required volume for 3 functionalizations: **100 µL** with a final concentration of **100 nM**.

Vial 1		Vial 2
<b>Adapter strand 1 - Ga</b> (400 nM)	Conjugated <b>Ligand strand</b> (500 nM)	<b>Adapter strand 2 - Ga - lfs</b> (200/250 nM)
25 µL	25 µL	50 µL

After incubation time, mix vial 1 and vial 2 to obtain 100 µL of ready-to-use DNA solution.

## Contact

**Dynamic Biosensors GmbH**

Perchtinger Str. 8/10  
81379 Munich  
Germany

**Dynamic Biosensors, Inc.**

300 Trade Center, Suite 1400  
Woburn, MA 01801  
USA

**Order Information**   [order@dynamic-biosensors.com](mailto:order@dynamic-biosensors.com)

**Technical Support**   [support@dynamic-biosensors.com](mailto:support@dynamic-biosensors.com)

**[www.dynamic-biosensors.com](http://www.dynamic-biosensors.com)**

Instruments and chips are engineered and manufactured in Germany.

©2024 Dynamic Biosensors GmbH | Dynamic Biosensors, Inc. All rights reserved.

[1] TE40: 10 mM Tris, 40 mM NaCl, 0.05 % Tween20, 50 µM EDTA, 50 µM EGTA

[2] If the protein is not stable in PE40 (TE40, HE40), please check buffer compatibility with the [switchSENSE®](#) compatibility sheet.