


Product Description

Product Code	ENZ-80-1-Y1-S
Suitable for	switchSENSE® DRX ^{yellow}
Storage	Store at 2-8 °C, dry in the dark
Layout	2 reference spots (nanolever A) and 4 probe spots (nanolever B)
Grade	Standard

Multi-purpose chips are versatile tools to analyze protein-protein, protein-small molecule, DNA/RNA-protein and DNA/RNA-small molecule interactions.

Biochip Design

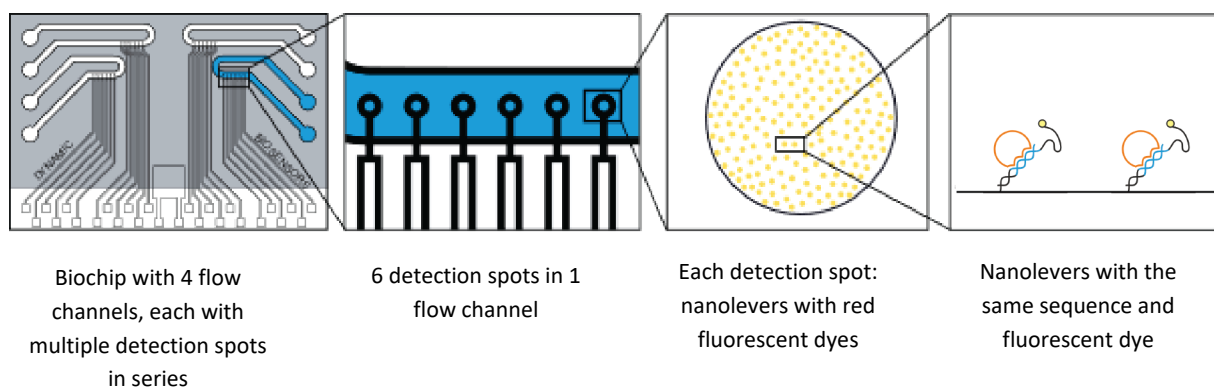

The detection spots are functionalized with one nanolever, NL-T80. The 3'-end of the DNA is fixed on the surface while the 5'-end carries a yellow fluorescent dye.

Biochip Layout | Flow channels with 6 single sequence (T) spots

With 6 single sequence spots, the electrode functionalization features replicates in favor of better measurement statistics.



For more information about the DNA sequences, please contact us at info@dynamic-biosensors.com.

DNA for the Biochip

Each Biochip comes with unmodified complementary DNA:

cNL-P48T80 (500 µL, c = 1 µM)

cNL-P80T80 (500 µL, c = 1 µM)

Application Areas

<i>Binding Kinetics*</i>	k_{on} , k_{off} , K_D
<i>Binding Affinity*</i>	K_D , titration curve, n (Hill coefficient)
<i>Protein Diameter*</i>	D_H (Hydrodynamic diameter)
<i>Conformational Change*</i>	ΔD_H (relative change of hydrodynamic diameter)
<i>Melting & Thermodynamics*</i>	T_M , ΔG , ΔH , ΔS
<i>Multimers & Aggregation*</i>	Monomer-dimer discrimination, aggregation
Nuclease & Polymerase Activity	k_{cat} , K_M , $T_{activate}$
<i>Bispecific Binders & Avidity*</i>	Binding affinity/avidity/kinetics

* for this application, other biochips (MPC or BIF) are recommended

Biochip Handling

The biochip is ready to use. For research use only.

Avoid touching or picking up the biochip with your hands as this may cause electrostatic discharge, which harms the bilayer. Instead, always use the vacuum tweezers provided with the instrument when handling the biochips.

After installing the biochip in the DRX instrument and before starting an experiment be sure that the selected channel has been passivated with passivation solution (Order No. SOL-PAS-1-5).

More information about biochip storage and handling can be found in the document “switchSENSE® Biochips Storage Information” – [DOWNLOAD HERE](#).

Grade I Standard

Standard grade chips are batch tested, which means selected chips undergo comprehensive quality control (switching dynamics, voltage-response of the DNA layer, and fluorescence levels). In addition, the fluorescence levels on each electrode are tested to be within specifications. Due to the production process not all detection spots fulfill our high quality standards. It is guaranteed that 20 or more detection spots are QC approved.

Compatible Functionalization Kits

All compatible functionalization kits for this biochip are listed on our website:

www.dynamic-biosensors.com/reagents/

Contact

Dynamic Biosensors GmbH
Lochhamerstr. 15
82152 Martinsried/Planegg
Germany

Dynamic Biosensors Inc.
9705 Carroll Centre Road, Suite 100
San Diego, CA 92126
USA

Order Information

Phone: +49 89 89 74 544 0
order@dynamic-biosensors.com

Technical Support

Phone: +49 89 89 74 544 66
support@dynamic-biosensors.com