


10x Buffer TE40 pH 7.4

running buffer

Dynamic Biosensors GmbH & Inc.

BU-TE-40-10 v2.1

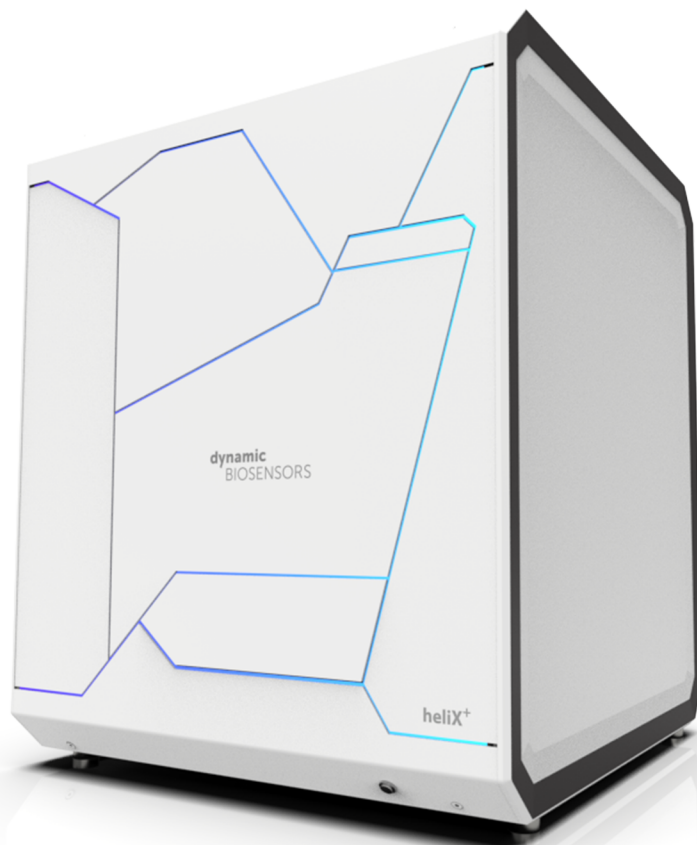

Product Description

Order Number: **BU-TE-40-10**

Table 1. Contents and Storage Information

Material	Composition	Amount	Storage
10x Buffer TE40 pH 7.4	100 mM Tris, 400 mM NaCl, 500 µM EDTA, 500 µM EGTA and 0.5 % Tween20; 0.2 µm sterile filtered	50 mL	2-8°C

For research use only.

This product has a limited shelf life, please see expiry date on label.

Preparation

Dilute the complete solution 10x Buffer TE40 pH 7.4 (50 mL) by mixing with 450 mL ultrapure water.

After dilution TE40 Buffer is ready for use (10 mM Tris, 40 mM NaCl, 50 µM EDTA, 50 µM EGTA and 0.05 % Tween20).

The diluted buffer should be stored at 2-8°C.

Contact

Dynamic Biosensors GmbH

Perchtinger Str. 8/10
81379 Munich
Germany

Dynamic Biosensors, Inc.

300 Trade Center, Suite 1400
Woburn, MA 01801
USA

Order Information order@dynamic-biosensors.com

Technical Support support@dynamic-biosensors.com

www.dynamic-biosensors.com

Instruments and chips are engineered and manufactured in Germany.

©2024 Dynamic Biosensors GmbH | Dynamic Biosensors, Inc. All rights reserved.